



Product Highlights

- ► Automatically records peak force and calculates CV and bilateral strength deficits
- ► PowerTrack II's thin profile improves stability during tests and "Axis Compensation" circuitry guarantees 99% accuracy for off center placement
- PowerTrack II's innovative pop-in "quick connect" system makes changing test pads faster and easier

Number of tests stored:

Tests bilaterally:

Tests unilaterally:

Number of repetitions per test:

Audible cues:

Load cell capacity:

Automatic coefficient of

variation calculations (CV):

Automatic bilateral deficit calculation:

Threshold settings:

Unit of measure (pounds default):

Test times:

Allows for specific test re-testing:

Print options:

Upgrades to Tracker Legacy computerized testing without modification:

Accessories:

1-40 Yes

Yes

1-4

Yes

125 lbs.

Yes

Yes

1, 3, 5, 10, or 20 lbs.

lb. or N.

Untimed


Yes

Yes

Yes

Protective case, flipchart, worksheets, two test pads and manual



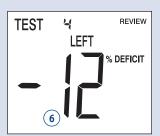

- 1 Test number
- 2 Repetition number
- 3 Active "mode" for testing or reviewing data
- 4 Side being tested
- 5 Real-time force on dynamometer shown in fractions of a pound or newtons
- 6 Bilateral comparison





Optional **DownLoader Software and Adapter** downloads data from your Commander Console into a computer text file for quick printing or easy inclusion into your documentation software.



Tracker™ MDE Software option allows you to document findings, track progress and quickly print reports for all Commander products.

Normative comparisons
CV calculations
AMA validity for spine ROM
Bilateral comparisons
Document sensory dermatomes, orthopedic tests, deep tendon reflexes
Score outcomes questionnaires
Comprehensive intake section
Optional AMA-licensed Impairment Module
More



(800) 985-8324 fax: (801) 478-0674 (801) 478-0680 www.jtechmedical.com